


SSCY_3X150 0.8

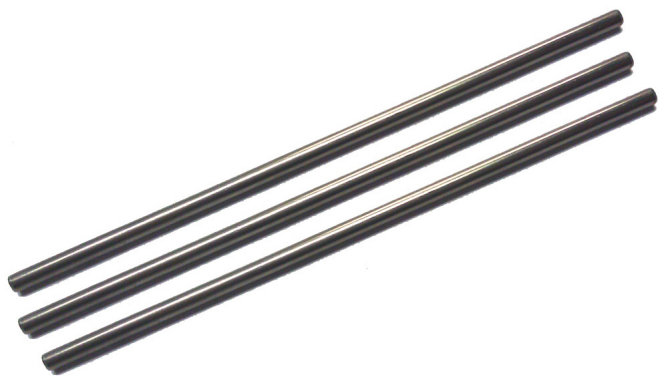

Stainless Steel cylinder D=0.8, Length=150 mm - 3 pcs

Rating: Not Rated Yet

Price

Salesprice with discount

23,00 €

3-5 Days

[Ask a question about this product](#)

Manufacturer [Linari Engineering srl](#)

Description

The drum collectors are made of high resistance hollow aluminum to reduce to a minimum both the weight and the moment of inertia so that they can be used at high speeds of rotation without triggering dangerous vibrations.

This collector and all the ones with length up to 150 mm can be used on all the rotation systems of Linari Engineering.

Every drum collector comes with bare, non-anodized and untreated aluminum outer surfaces to reduce any possible electrical resistance which would alter the electrical field making the electrospinning process non-reproducible.

Dimensions

Outer diameter	0.83 mm
Inner diameter	0.51 mm
Lenght	150 mm

Packaging

3 pcs/conf

Units in box: 1